

CHAPTERS

Coherent Sources

Incoherent Sources

Quantum Electronics

Drivers/Mounts

Accessories

SECTIONS

Laser Glasses

Laser Safety Signs

Safety and Blackout

Beam Blocks

Alignment Tools

ESD Protection

For current pricing,  
please see our website.

# Laser Safety Signs

The LSL10 laser safety light box provides a clear and concise indication that a laser system is in use, illuminating the appropriate user-chosen warning message and laser classification. Designed to meet ANSI Z136.1, the ANSI Standard for Safe Use of Lasers, the box operates from 110 VAC or 220 VAC (LSL10-EC) and has safety interlock features that prevent the use of interlock-equipped laser systems unless the safety light is turned on.

## Features

- Provides a Clear and Concise Indication that a Laser is in Use
- Four Signs Available
- Dimensions: 10" x 12" x 3" (254 mm x 305 mm x 76 mm)

Measuring 10" x 12" (254 mm x 305 mm), these illuminated signs have large, legible text. The LSL10 series of signs are useful tools for any research or manufacturing facility.


The LSL10 accepts any one of the four signs available for covering Laser Classes 2 through 4. Each sign comes with a set of alpha-numeric labels enabling each user the ability to customize the sign for a specific wavelength and output power as required by ANSI Z136.1 and other pertinent laser safety specifications.

These laser safety signs can be displayed on desktops or benchtops, or they can be mounted to any wall. A convenient line cord latch prevents the AC line cord from disengaging from the unit.


Complies with ANSI-Z136.1




LSL10 Shown with LSL10-3A  
10" x 12" x 3" (254 mm x 305 mm x 76 mm)

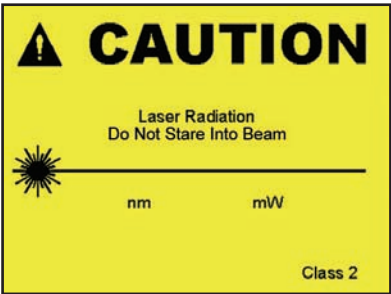


LSL10-2

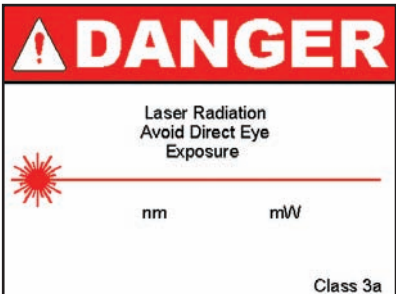

LSL10-3A

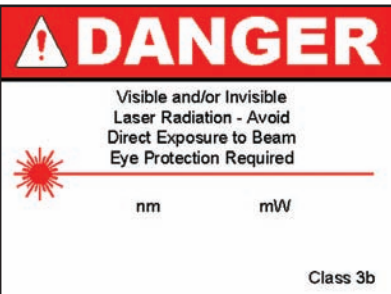

LSL10-3B

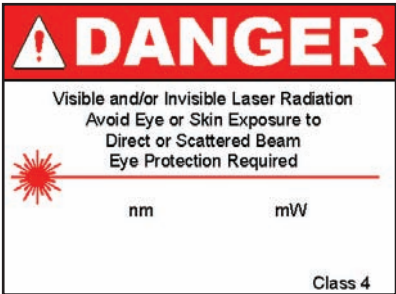

LSL10-4

ITEM #	\$	£	€	RMB	DESCRIPTION
LSL10	\$ 265.00	£ 190.80	€ 230,55	¥ 2,112.05	Lighted Laser Safety Box, 110 VAC Input*
LSL10-EC	\$ 265.00	£ 190.80	€ 230,55	¥ 2,112.05	Lighted Laser Safety Box, 220 VAC Input*
LSL10-2	\$ 48.00	£ 34.56	€ 41,76	¥ 382.56	Class 2 Laser Safety Sign
LSL10-3A	\$ 48.00	£ 34.56	€ 41,76	¥ 382.56	Class 3A Laser Safety Sign
LSL10-3B	\$ 48.00	£ 34.56	€ 41,76	¥ 382.56	Class 3B Laser Safety Sign
LSL10-4	\$ 48.00	£ 34.56	€ 41,76	¥ 382.56	Class 4 Laser Safety Sign

\*Signs not included