


## NR Series SM1-Compatible Adapter Plate



NR360SP8

This adapter plate bolts to the top of the NR360S NanoRotator™ Stage. An SM1-threaded through hole allows Ø1" optics and Thorlabs' SM1 family of products to be mounted to the stage. Thorlabs' standard 30 mm cage assemblies can be attached via 4-40 taps.

- SM1-Threaded Through Hole
- Four Each of 4-40 , 8-32 and, 1/4"-20 Taps
- #8 (M4) Counterbores for Securing the Plate to the Stage

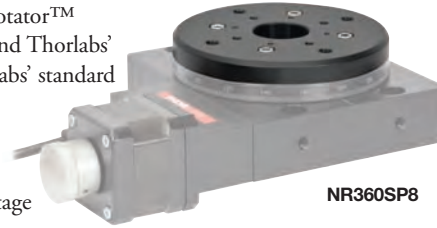

NR360SP8

ITEM #	\$	£	€	RMB	DESCRIPTION
NR360SP8	\$ 48.00	£ 34.56	€ 41,76	¥ 382.56	SM1-Compatible Adapter Plate for the NR360S Stage

## NR Series Threaded Adapter Plate



NR360SP9

This adapter plate is used to provide a number of different mounting options when attaching devices to our NR360S rotation stage. The carriage has 17 1/4"-20 (M6 x 1.0) through holes. The four counterbores for #8 (M4) cap screws allow the plate to be attached to the NR360S (NR360S/M) Stage.

**NEW**  
product



NR360SP9

ITEM #	METRIC ITEM #	\$	£	€	RMB	DESCRIPTION
NR360SP9	NR360SP9/M	\$ 88.00	£ 63.36	€ 76,56	¥ 701.36	Threaded Adapter Plate for the NR360S Stage

## NR Series Grooved Adapter Plate



Adapter Plate Shown with Waveguide and Fiber Accessories (See Pages 578 – 582)

This adapter plate provides an easy means to integrate all of the miscellaneous accessories from our fiber launch product line, including optical mounts, tip-tilt mounts, and fiber holders. Please refer to page 572 for the complete selection of products.



NR360SP4

ITEM #	METRIC ITEM #	\$	£	€	RMB	DESCRIPTION
NR360SP4	NR360SP4/M	\$ 112.00	£ 80.64	€ 97,44	¥ 892.64	Grooved Adapter Plate for the NR360S Stage

## Vertical Mounting Bracket

The NR360SP2 Vertical Mounting Bracket allows the NR360S NanoRotator™ Stage to be mounted in a vertical orientation. The bracket can be attached to an optical table or breadboard using the seven 1/4"-20 (M6 x 1.0) mounting holes.

The photograph to the right shows the NR360SP2 with a NanoRotator NR360S.



NR360SP2



ITEM #	METRIC ITEM #	\$	£	€	RMB	DESCRIPTION
NR360SP2	NR360SP2/M	\$226.00	£162.72	€ 196,62	¥ 1,801.22	Vertical Mounting Bracket

## NR Series Adapter Plates

The NR360SP5 allows the NR360S NanoRotator™ to be mounted directly to any LNR or NRT stage. The LNR stage is found on pages 514 – 516, while the NRT stage is found on pages 518 – 519.

NR360SP5  
Bottom Mounting  
Adapter



ITEM #	METRIC ITEM #	\$	£	€	RMB	DESCRIPTION
NR360SP5	NR360SP5/M	\$ 89.00	£ 64.08	€ 77,43	¥ 709.33	Mounting Adapter Plate