

Ø1/2" Post Holders

Each Ø1/2" Post Holder has a wide, square relief cut that provides a highly stable two-line contact for the mating post. A spring-loaded plunger within the brass thumbscrew holds the steel post firmly within the post holder prior to final tightening. This allows easy, one-handed height adjustment of post-mounted devices.


ITEM#	LENGTH	METRIC ITEM#	LENGTH	\$	£	€	RMB
PH1	1"	PH1/M	25.4 mm	\$ 7.03	£ 4.85	€ 6.30	¥ 59.40
PH1.5	1.5"	PH1.5/M	38.1 mm	\$ 7.22	£ 5.00	€ 6.50	¥ 61.00
PH2	2"	PH2/M	50.8 mm	\$ 7.70	£ 5.35	€ 6.90	¥ 65.10
PH3	3"	PH3/M	76.2 mm	\$ 8.27	£ 5.75	€ 7.40	¥ 69.90
PH4	4"	PH4/M	101.6 mm	\$ 9.17	£ 6.35	€ 8.20	¥ 77.50
PH6	6"	PH6/M	152.4 mm	\$ 12.64	£ 8.75	€ 11.30	¥ 106.80



PH-Series
Post Holders

**NEW
design**

Please refer to our website for complete models and drawings.

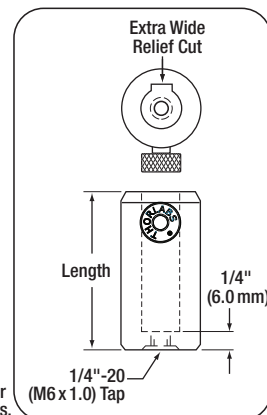


Ø1/2" Translating Post Holders

Our Translating Post Holders provide lockable, non-rotating height adjustment with convenient one-handed operation. A brass thumbscrew with a spring-loaded plunger provides holding force while the wide, square relief cut within the Ø1/2" post bore offers a stable two-line contact. A knurled adjustment ring actuates the non-rotating height adjustment, offering 0.025" (1.0 mm) of fine translational adjustment per revolution.

Once a desired position is achieved, the entire mechanism can be locked in place with the lower thumbscrew.

Please refer to our website for complete models and drawings.

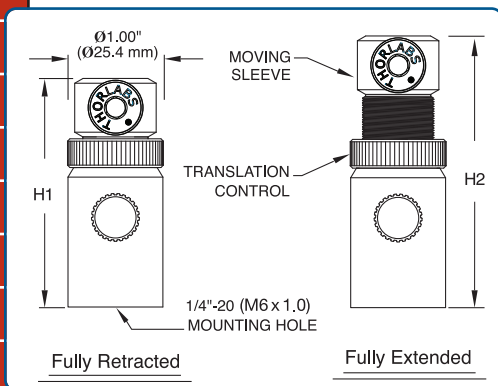

PH2T

PH3T

	H1	H2
PH2T	2.38" (60.5 mm)	2.83" (71.9 mm)
PH3T	2.89" (73.4 mm)	3.53" (89.7 mm)

ITEM#	METRIC ITEM#	\$	£	€	RMB	SIZE
PH2T	PH2T/M	\$ 46.00	£ 31.90	€ 40.90	¥ 388.50	2" (50.8 mm) Translating Post Holder
PH3T	PH3T/M	\$ 47.90	£ 33.30	€ 42.60	¥ 404.50	3" (76.2 mm) Translating Post Holder

Locking Thumbscrew for Post Holders




TS25B031

This brass thumbscrew with hex broach offers extra holding power while protecting the steel post from marring. The TS25B031 is compatible with both imperial and metric post holders.

ITEM#	\$*	£*	€*	RMB*
TS25B031	\$ 12.40	£ 8.60	€ 11.10	¥ 104.80

*Priced and Packaged 10 Per Box

Spring-Loaded Thumbscrew

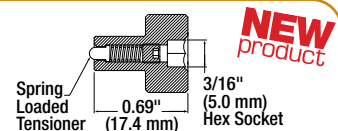

TS25H
Spring-Loaded
Tip

The TS25H is designed to firmly hold steel posts within post holders prior to final tightening, allowing one-handed positioning of post-mounted devices.

ITEM#*	\$	£	€	RMB
TS25H	\$ 2.10	£ 1.45	€ 1.90	¥ 17.80


*Compatible with Imperial and Metric Post Holders

Swivel Base Adapter



UPHA

Our 360° Swivel Base Adapter fits all standard PH-Series Post Holders with 1/4" (M6) through holes. The swivel design is perfect for applications that require a large number of components to be packed into a small area. The swivel base decreases the time required to assemble an optical system by allowing a multitude of mounting positions. Magnets located underneath the base facilitate temporary positioning prior to fastening.

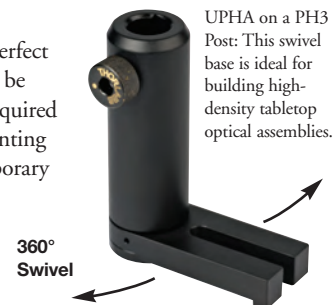

ITEM#	\$	£	€	RMB	DESCRIPTION
UPHA	\$ 19.90	£ 13.80	€ 17.70	¥ 168.10	Universal Post Holder Adapter